
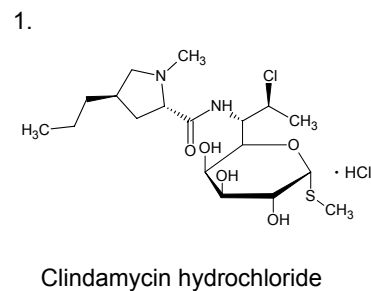


クリンダマイシン塩酸塩（第十五改正日本薬局方第二追補記載条件）

Clindamycin hydrochloride (The Japanese Pharmacopoeia 15th Edition Supplement II) F100416B

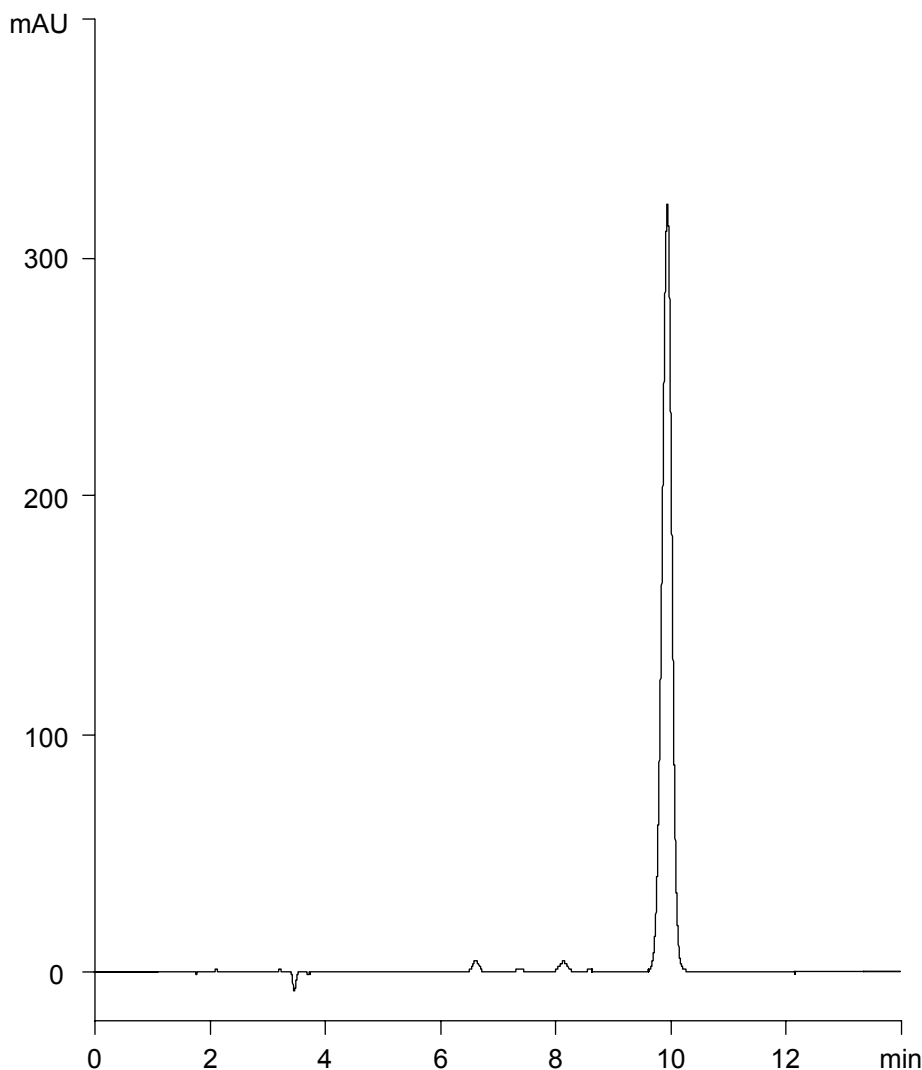

Column : YMC-Triart C18 (5 μ m, 12 nm)
250 X 4.6 mm I.D.
Eluent : 50 mM KH_2PO_4 (pH 7.5 adjusted by 8 M KOH)/acetonitrile (55/45)
Flow rate : 1.1 mL/min
Temperature : 25°C
Detection : UV at 210 nm
Injection : 20 μ L (1.0 mg/mL)

(第十五改正日本薬局方第二追補記載条件；定量法)