


カンナビノイド
Cannabinoids

R170328A

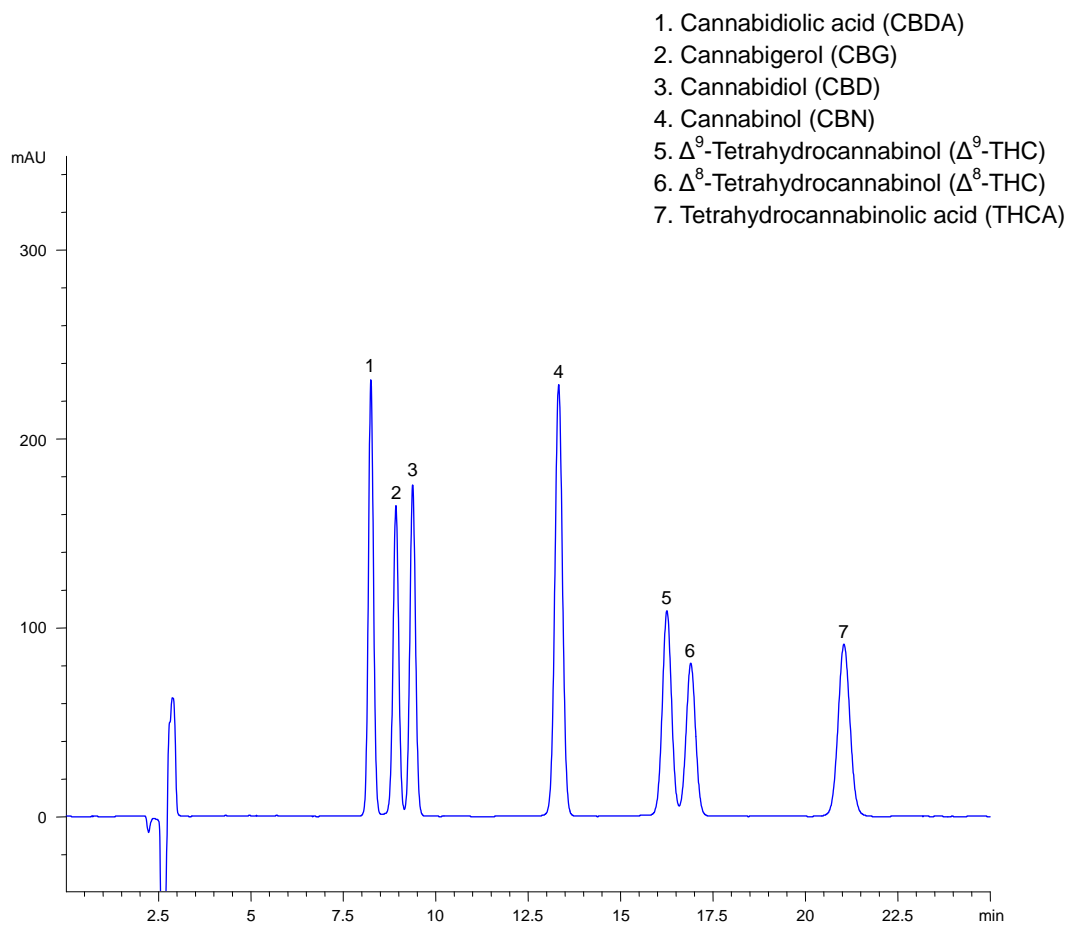


Courtesy of YMC America, Inc.

Column	: YMC-Triart C18 (5 μ m, 12 nm) 250 X 4.6 mmI.D.
Eluent	: acetonitrile/water/acetic acid (79/21/0.1)
Flow rate	: 1.0 mL/min
Temperature	: 35°C
Detection	: UV at 220 nm
Injection	: 25 μ L (0.05 mg/mL)

カンナビノイド
Cannabinoids

R170328B

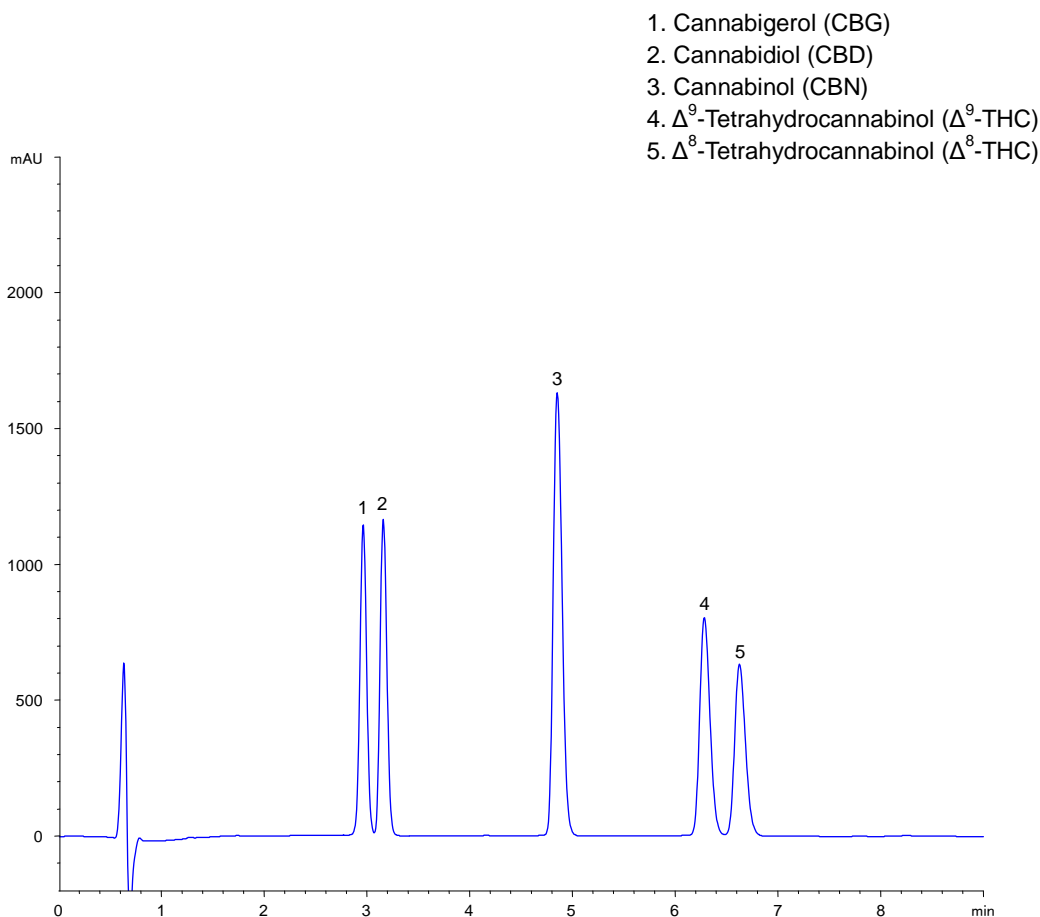


Courtesy of YMC America, Inc.

Column	: Meteoric Core C18 (2.7 μ m, 8 nm) 100 X 4.6 mmI.D.
Eluent	: acetonitrile/water/formic acid (73/27/0.1)
Flow rate	: 1.25 mL/min
Temperature	: 35°C
Detection	: UV at 220 nm
Injection	: 25 μ L (0.2 mg/mL)

カンナビノイド
Cannabinoids

R170328C

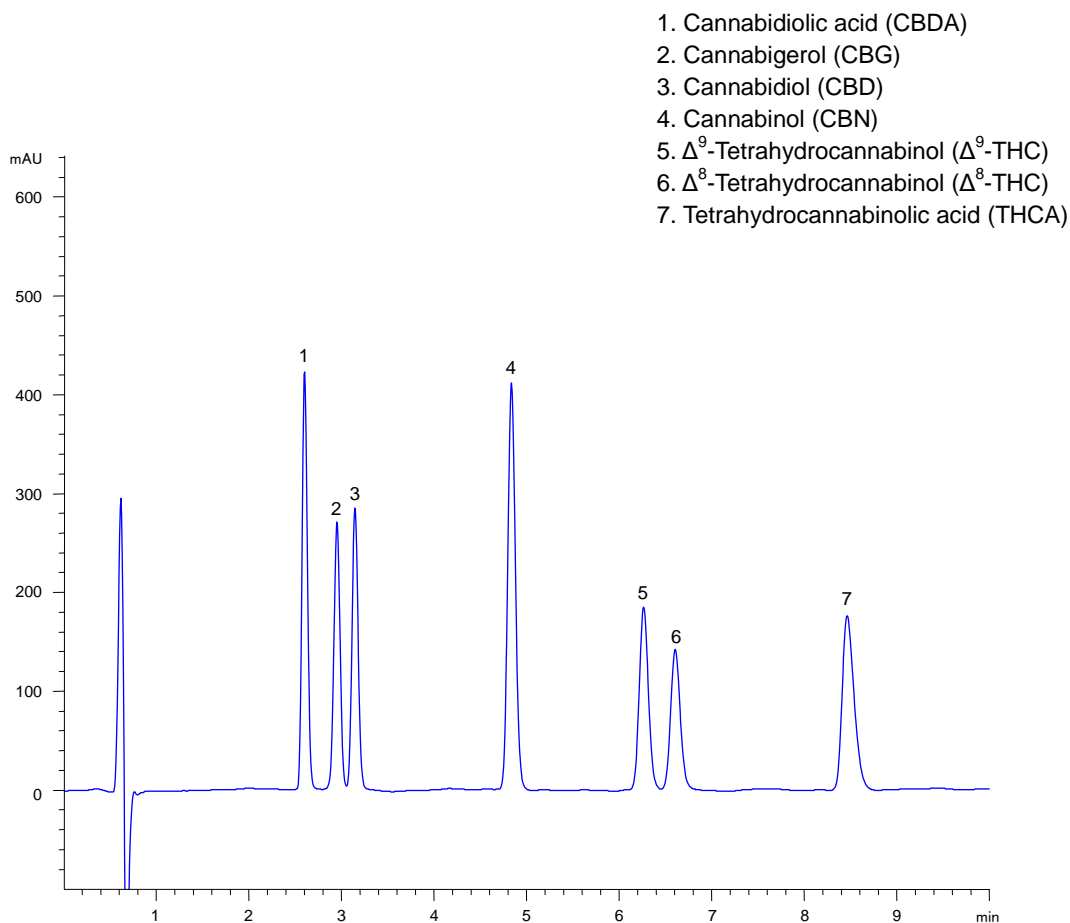

Courtesy of YMC America, Inc.

Column	: Meteoric Core C18 (2.7 μ m, 8 nm) 100 X 4.6 mm I.D.
Eluent	: acetonitrile/water/formic acid (73/27/0.1)
Flow rate	: 1.25 mL/min
Temperature	: 35°C
Detection	: UV at 220 nm
Injection	: 25 μ L (0.05 mg/mL)