


抗体薬物複合体 (ブレンツキシマブ ベドチン)
Antibody-Drug Conjugates (Brentuximab vedotin)

C240226B

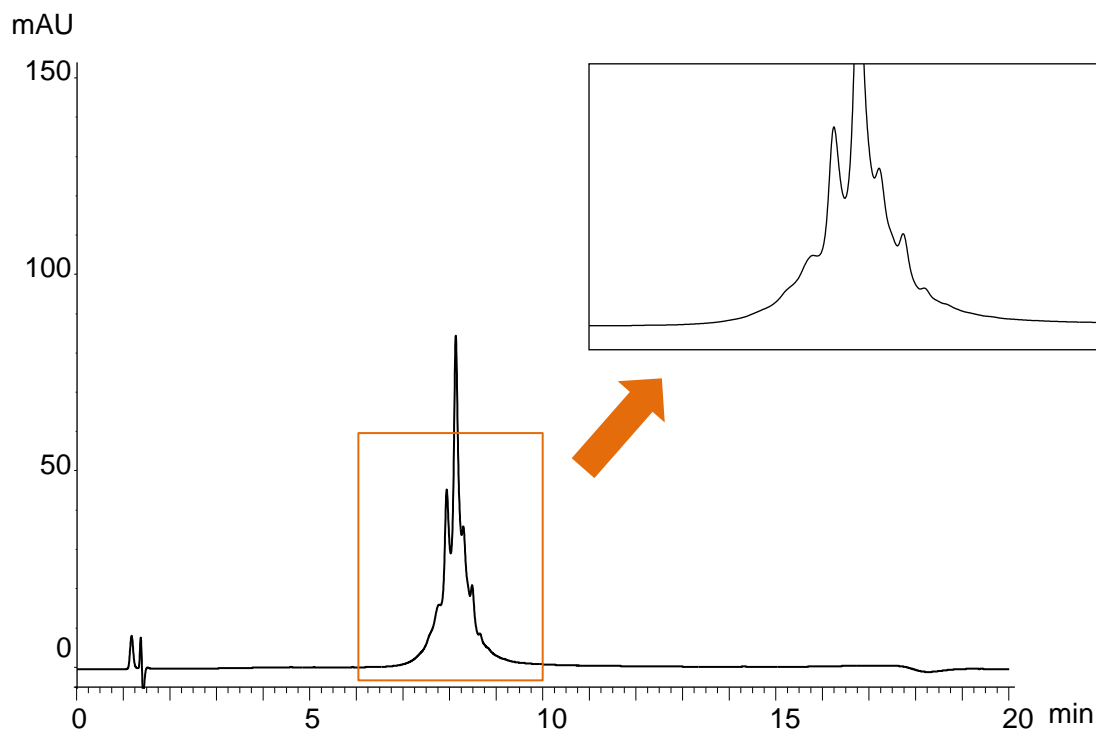

Column : Accura BioPro IEX SF (3 μ m)
100 X 4.6 mm I.D.

Eluent : A) 10 mM MES*-NaOH (pH 6.0)
B) 10 mM MES*-NaOH (pH 6.0) containing 1.0 M NaCl
0–20%B (0–15 min)

Flow rate : 0.5 mL/min

Temperature : 25°C

Detection : UV at 280 nm

Injection : 10 μ L (1 mg/mL)

Sample : Brentuximab vedotin for injection
Courtesy of Prof. S. Manabe, Hoshi University, Japan/Tohoku University, Japan

*2-(N-morpholino) ethanesulfonic acid