

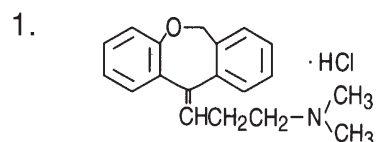


三環系抗うつ薬
Tricyclic antidepressants


D960724A

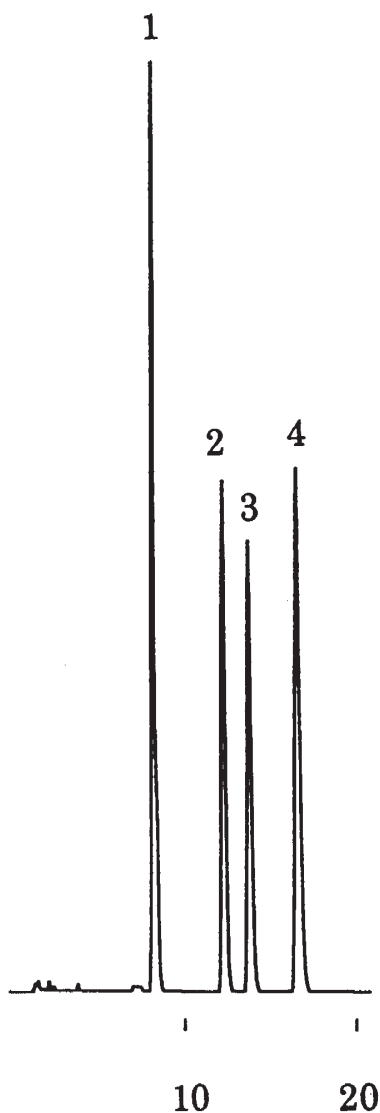


Doxepin hydrochloride

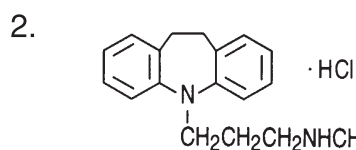

Desipramine hydrochloride

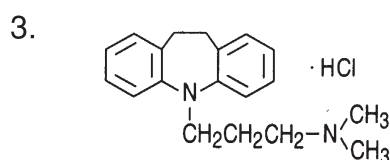

Imipramine hydrochloride

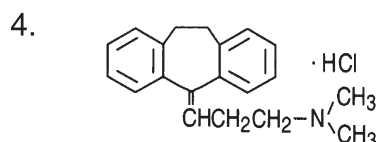

Amitriptyline hydrochloride

Column : YMC-Pack Pro C18 (5 μm, 12nm)
150 × 4.6 mm I.D.
Eluent : acetonitrile/50mM NH₄H₂PO₄ (30/70)
Flow rate : 1.0 mL/min
Temperature : 37 °C
Detection : UV at 215 nm, 0.13 AUFS
Injection : 2 μL (0.125 mg/mL)