



抗体医薬品
Antibody drugs

U240403A

Bevacizumab

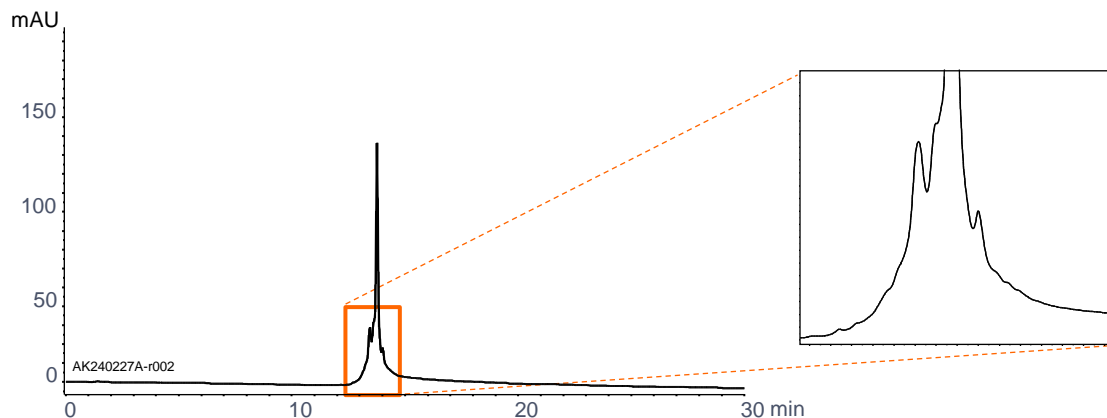

Trastuzumab

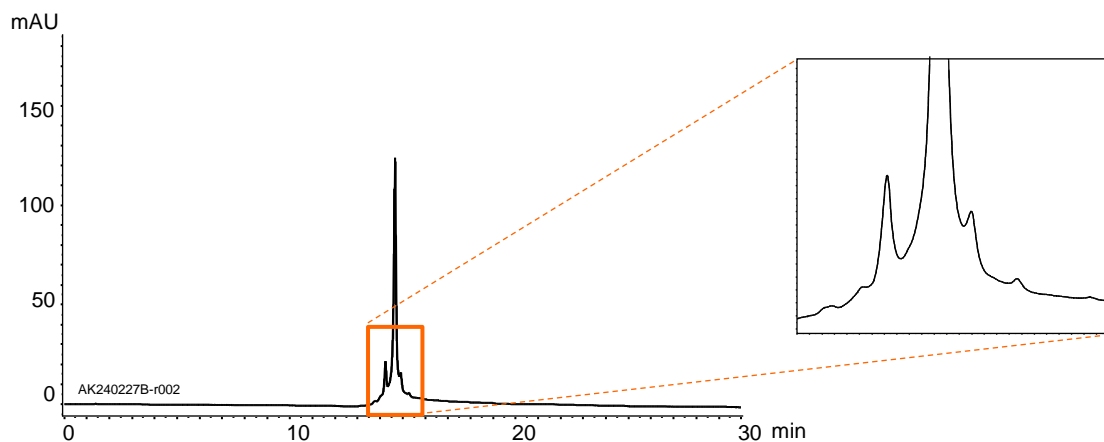

Column : Accura BioPro IEX SF (3 μ m)
100 X 2.1 mm I.D.

Eluent : A) 10 mM MES*-NaOH (pH 6.6)
B) 10 mM MES*-NaOH (pH 6.6) containing 1.0 M NaCl
0-20%B (0-15 min)

Flow rate : 0.1 mL/min

Temperature : 25°C

Detection : UV at 280 nm


Injection : 2 μ L (1 mg/mL)

*2-(*N*-morpholino) ethanesulfonic acid


抗体医薬品
Antibody drugs

U240403B

Bevacizumab

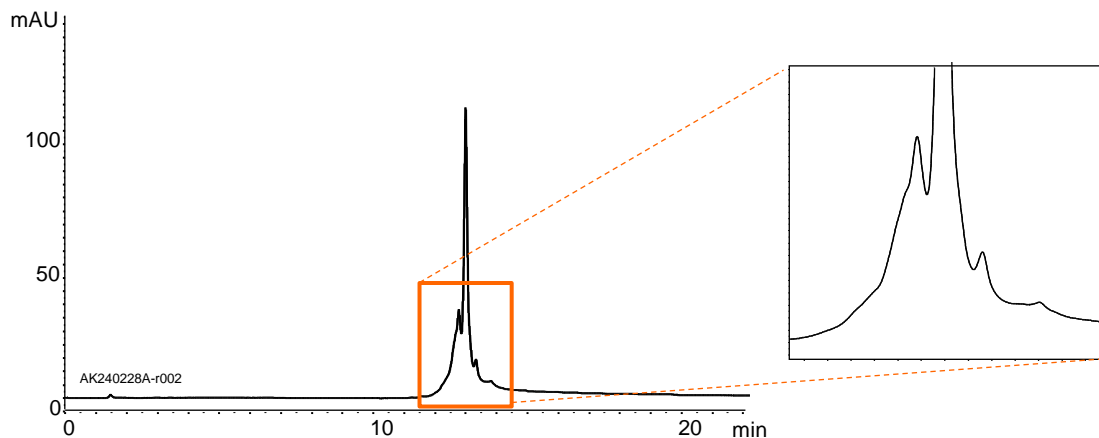

Trastuzumab

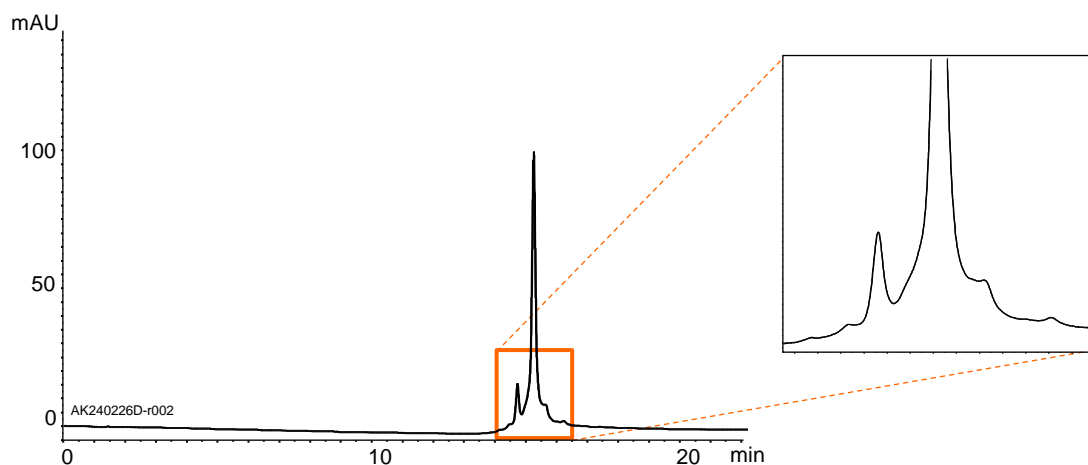

Column	: Accura BioPro IEX SF (3 μ m) 100 X 2.1 mm I.D.
Eluent	: A) 20 mM $\text{CH}_3\text{COONH}_4$ - CH_3COOH (pH 5.6) B) 140 mM $\text{CH}_3\text{COONH}_4$ -10 mM NH_4HCO_3 (pH 7.4) 20%B (0–2 min), 20–100 %B (2–18 min), 100 %B (18–22 min)
Flow rate	: 0.1 mL/min
Temperature	: 25°C
Detection	: UV at 280 nm
Injection	: 2 μ L (1 mg/mL)